
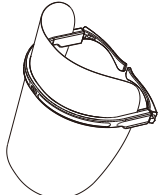
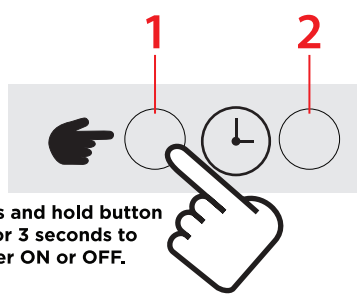
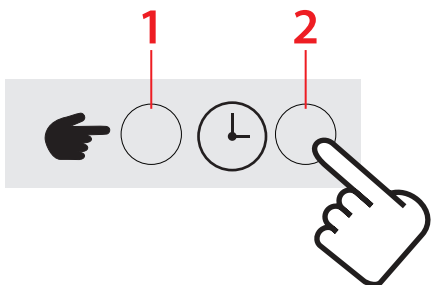


# Lomi User Guide

## (Front)

<h3>Compliance Information</h3> <p><b>FCC Information</b> This device complies with Part 15 of the FCC Rules.</p> <p><b>Caution:</b> Any changes or modifications not expressly approved could void the user's authority to operate the equipment.</p> <p><b>Note:</b> This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:</p> <ul style="list-style-type: none"> <li>• Reorient or relocate the receiving antenna.</li> <li>• Increase the separation between the equipment and receiver.</li> <li>• Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.</li> <li>• Consult the dealer or an experienced radio/TV technician for help.</li> </ul> <p>This device complies with Part 15 of the FCC Rules, Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.</p> <p><b>Radiation Exposure Information</b> This device has been evaluated to meet general RF exposure requirements. The device can be used in portable exposure conditions without restriction.</p> <p><b>WARNING Li-Ion Battery Inside</b> This product is fitted with a Li-Ion battery. Do not damage, open, or dismantle the battery and do not use it in damp and/or corrosive conditions. Use only with compatible chargers. Never dispose of batteries in a fire, and never expose them to high temperatures. Do not expose the product to temperatures exceeding 60°C (140°F).</p> <p><b>Rating Label Definitions</b></p> <ul style="list-style-type: none"> <li> Direct current symbol.</li> <li> This symbol signifies that the product must not be discarded as household waste, and should be delivered to an appropriate collection facility for recycling. Proper disposal and recycling helps protect natural resources, human health and the environment.</li> <li> This symbol signifies and declares that this equipment complies with battery charging which complies with California CEC and energy efficiency.</li> <li> This symbol represents FCC Declaration of Conformity and compliance.</li> </ul>	<h3>Warranty Information</h3> <p><b>90-Day Limited Warranty</b> Territory: United States/Canada</p> <p><b>LIMITED WARRANTY TO ORIGINAL CONSUMER</b> This product as supplied and distributed new by an authorized retailer is warranted by Southern Telecom, Inc. to the original consumer purchaser against defects in materials and workmanship ("Warranty") as follows:</p> <p><b>To obtain warranty service:</b></p> <ul style="list-style-type: none"> <li>• Visit our customer support portal: <a href="http://www.customersupport123.com">www.customersupport123.com</a></li> <li>• Select the brand and model of your device</li> <li>• Select "Service Request" and fill out the form to begin your inquiry</li> <li>• Exchanges require dated proof of purchase from an authorized retailer</li> </ul> <p>Your product will be repaired or replaced, at our option, for the same or similar model of equal value if examination by the service center determines this product is defective. Products received damaged as a result of shipping will require you to file a claim with the carrier.</p> <p><b>Warranty service not provided</b> This warranty does not cover damage resulting from accident, misuse, abuse, improper installation or operation, lack of reasonable care, unauthorized modification including software modifications such as the installation of custom firmware. This warranty is voided in the event any unauthorized person opens, alters or repairs this product. All products being returned to the authorized service center for repair must be suitably packaged.</p> <p><b>NO WARRANTIES, WHETHER EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OTHER THAN THOSE EXPRESSLY DESCRIBED ABOVE SHALL APPLY. DISTRIBUTOR FURTHER DISCLAIMS ALL WARRANTIES AFTER THE EXPRESS WARRANTY PERIOD STATED ABOVE. NO OTHER EXPRESS WARRANTY OR GUARANTY GIVEN BY ANY PERSON, FIRM OR ENTITY WITH RESPECT TO THE PRODUCT SHALL BE BINDING ON DISTRIBUTOR, REPAIR, REPLACEMENT, OR REFUND OF THE ORIGINAL PURCHASE PRICE - AT DISTRIBUTOR'S SOLE DISCRETION - ARE THE EXCLUSIVE REMEDIES OF THE CONSUMER. IN NO EVENT WILL DISTRIBUTOR OR ITS MANUFACTURERS BE LIABLE FOR ANY INCIDENTAL, DIRECT, INDIRECT, SPECIAL, PUNITIVE OR CONSEQUENTIAL DAMAGES (SUCH AS, BUT NOT LIMITED TO, DAMAGES FOR LOSS OF PROFITS, BUSINESS, SAVINGS, DATA OR RECORDS) CAUSED BY THE USE, MISUSE OR INABILITY TO USE THE PRODUCT. EXCEPT AS STATED HEREIN, NO OTHER WARRANTIES SHALL APPLY. NOTWITHSTANDING THE FOREGOING, CONSUMER'S RECOVERY AGAINST DISTRIBUTOR SHALL NOT EXCEED THE PURCHASE PRICE OF THE PRODUCT SOLD BY DISTRIBUTOR. THIS WARRANTY SHALL NOT EXTEND TO ANYONE OTHER THAN THE ORIGINAL CONSUMER WHO PURCHASED THE PRODUCT AND IS NOT TRANSFERABLE.</b></p> <p>Some countries, states or provinces do not allow the exclusion or limitation of incidental or consequential damages or allow limitations on warranties, so limitation or exclusions may not apply to you. This warranty gives you specific legal rights, and you may have other rights that vary from state to state or province to province. Contact your authorized retailer to determine if another warranty applies.</p> <p><b>All product names, trademarks, and registered trademarks are property of their respective owners. All company, product, and service names used in this document are for identification purposes only. Use of these names, trademarks, and brands does not imply endorsement.</b></p> <p>Printed in China Designed in the U.S. v1, 12/2020</p>	 <p><b>LOMI</b> Light Therapy Facial Mask</p> <p>User Guide LOMB1103</p>
---	---	--

<h3>Before You Begin...</h3> <p><b>READ CAREFULLY AND KEEP THIS MANUAL</b></p> <p><b>NOT Intended for Medical Use</b> This product is designed as a non-professional device for skin care. <b>It is not intended to be a substitute for medical attention.</b> If your ailments persist or worsen after use, contact your doctor immediately.</p> <p><b>DO NOT USE IF:</b></p> <ul style="list-style-type: none"> <li>• You are pregnant.</li> <li>• You have any known thyroid issue.</li> </ul> <p><b>Please consult your doctor prior to using this product, if you have any concerns regarding your health.</b></p> <p><b>Caution</b></p> <ul style="list-style-type: none"> <li>• THIS IS NOT A TOY, <b>please keep away from children.</b></li> <li>• Keep away from open flames, water, moisture, and extreme heat or cold.</li> <li>• Supervision is recommended when this appliance is used by individuals with certain disabilities.</li> <li>• Use this appliance only for its intended use as described in this manual. DO NOT use attachments not recommended by the manufacturer.</li> <li>• Never operate this appliance if it is not working properly or if it has been dropped or damaged.</li> <li>• Keep the appliance away from heated surfaces.</li> <li>• <b>Keep away from the neck, Light should not be shined on the thyroid.</b></li> <li>• <b>Avoid shining light in the eyes. This can cause discomfort or injury.</b></li> <li>• <b>Do not use on sunburned, burned or damaged skin.</b></li> <li>• <b>Do not use on skin after applying medications, ointments or any other skin treatment that can make the skin more sensitive.</b></li> <li>• Do not open this product or attempt to repair the unit yourself.</li> </ul> <p><b>Questions &amp; Support</b> If you have any questions or problems concerning this product, please visit our website at <a href="http://customersupport123.com">customersupport123.com</a> and click on the <b>Brand</b> of your product.</p> <p><b>In the Box</b></p> <ul style="list-style-type: none"> <li>• Light Therapy Facial Mask (x1)</li> <li>• USB Power Cable (x1)</li> <li>• Instruction Manual (x1)</li> </ul> 	<h3>Getting Started</h3> <p><b>Charging the Battery</b> Plug the small connector of the included charging cable into the power port on the underside of the arm of the facial mask. Plug the larger connector into a USB power source, or into a <b>5V USB adapter (not included)</b> and plug into a wall socket.</p> <p>The LED indicator will <b>turn RED</b> while charging and turn OFF when fully charged.</p> <p><b>Powering ON/OFF</b> <b>Charge the facial mask fully before use.</b></p> <ul style="list-style-type: none"> <li>• Press and hold button (1) for 3 seconds to power ON and OFF.</li> </ul>  <p><b>Press and hold button (1) for 3 seconds to power ON or OFF.</b></p> <p><b>Using the Facial Mask</b></p> <ol style="list-style-type: none"> <li>1. Power ON. The lights will cycle through different colors every 5 minutes.</li> <li>2. Place the mask carefully on your clean, dry face, avoiding shining light into your eyes.</li> <li>3. Allow the mask to cycle through the different light colors on its own, or <b>press button (1) to switch between the 3 light colors.</b> <ul style="list-style-type: none"> <li>• <b>Red Light (620-750 nm Wavelength):</b> Lightens dark spots &amp; locks in moisture</li> <li>• <b>Yellow Light (590-620 nm Wavelength):</b> Increases elasticity &amp; restores balance</li> <li>• <b>Blue Light (476-495 nm Wavelength):</b> Treats acne &amp; rejuvenates skin</li> </ul> </li> </ol>	<p><b>Using the Timer</b> Press button (2) to turn the mask's timer on. Each press of button (2) will increase the timer by 5 minutes, with a maximum of 20 minutes. It is not recommended to use the facial mask for longer than 20 minutes.</p>  <p><b>Cleaning the Mask</b> The facial mask can be cleaned by wiping down the surfaces with alcohol pads or wipes. <b>Never submerge the mask in liquid,</b> or allow liquid to come into contact with the charging port, buttons, or charging cable.</p> <p><b>Storage</b> When it won't be used for a long period of time, the mask should be stored in a cool, dry location, away from sources of water, heat or dust. It is best stored in an airtight bag or container (<b>not included</b>).</p>
--	--	--

## (Back)