




LOMI

HANDHELD PORE VACUUM with 6 Head Attachments




User Guide
L0MB4004

Location of Controls

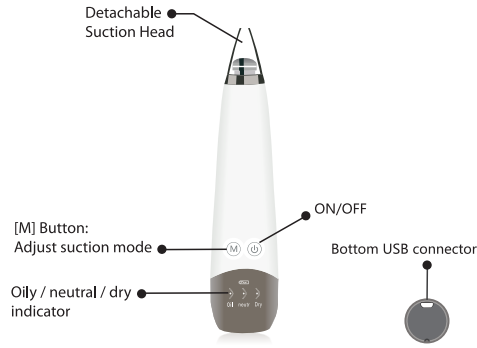

In the Box

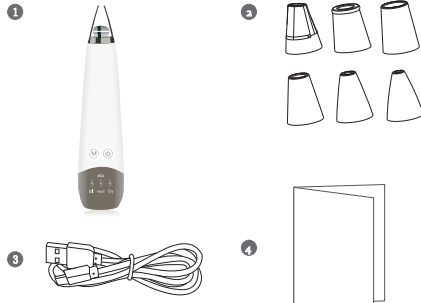

- 1 Handheld Pore Vacuum
- 2 Six suction heads
- 3 USB charging cable x1
- 4 Product manual x1

Introduction to Suction Heads

Six replaceable suction heads



Oval-shaped suction head

It is mainly used to remove fine lines. It needs to be quickly pulled away after sucking the skin for many times. It can promote the blood circulation of the skin and increase the elasticity of the skin.




Microcrystalline suction head

Fine microcrystalline head can gently remove rough old skin and aged toughness of the skin, remove dirt and debris, and also promote blood circulation and facilitate metabolism.



Large opening suction head

Strong suction. Do not use in one place for a long time. Move from bottom to top, moving slowly from the chin upward.

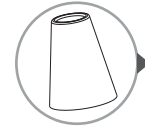

Medium opening suction head

Strong suction for blackheads or too much face oil. Do not use in one place for a long time. Move from bottom to top, moving slowly from the chin upward.



Small opening suction head

Light suction, for blackheads. For thin and tender skin or sensitive muscle, use this light suction head.



Mini opening suction head

Light suction, for sensitive skin.



Product on date	Executive standard: GB4706.1-2005 GB4706.15-2008
Quality inspector	Charging time: about 90 minutes
Certificate of approval	Suction : Up to 53kpa
	Rated voltage : 3.7V
	Input power : 2.5W
	Product model : MP-503
	Name: blackhead cleaner

Product specifications and electrical parameters

On low power, the battery LED indicator will flash RED. Please recharge when fully charged.
The battery LED indicator will flash RED while charging and turn OFF when fully charged.
adapter (not included) and plug into a wall outlet.
charging port on the pore vacuum. Plug the large connector into a USB power source, such as the USB port on a computer, or plug into a **5V USB**
Plug the small connector of the included charging cable into the

Charging the pore vacuum

Using on fine lines

1. Install the oval-shaped suction head.
2. Slowly use on the corners of the eyes, nose, lips and other parts. Adjust the suction strength as required. Pull up on the fine line repeatedly. This repeated action can penetrate skin deeply to promote firmness.

How to clean pore dirt

1. Install the large opening suction head.
2. Align on the suction head hole with the pore. Press slowly on the pore, and repeat several times.

Exfoliating the skin

1. Clean and dry the face.
2. Press on the microcrystalline or large opening suction head onto the pore vacuum.
3. At the bottom of the face, move from the lower jaw to the bottom. 4. On the cheeks, move from inside to outside.
5. On the forehead, start from the middle and move to the sides.
6. On the T-area, move from the bottom to the top.
7. This product is not suitable for people with physical, sensory or mental disability or lack of experience and knowledge including not allowed to be used by children, for example, the use must be supervised by the guardian from the side, and the use of this product be able to guide and take responsibility.
8. This product is not suitable for cosmetic surgery, nose augmentation, etc. - cut the parts that have undergone skin beauty surgery, pregnant women, it is not recommended for patients with heart disease.

Product usage

Please charge the pore vacuum fully before using.

Matters needing attention

1. Please dry the water or other liquid on your hand before use to avoid water or other liquid from the back of the body. The exhaust air hole enters the machine body, thus damaging the machine body.
2. Please do not use scrub or other harsh products before and after use.
3. With micro crystal head, dry skin - once a week, oily mixed skin every week can be used for 2 times, each time can not exceed 5 minutes, depending on the skin condition.
4. Please do not press the suction head on the skin, because the instrument will produce strong lifting force.
5. Please do not use the same part continuously. It may cause red and purple skin. If this happens, please suspend use.
6. After using, the skin will appear some irritation or redness due to the acceleration of internal circulation. This is a normal phenomenon. It will subside automatically about 30 minutes after use. Please keep skin warm.
7. This product is not suitable for people with physical, sensory or mental disability or lack of experience and knowledge including not allowed to be used by children, for example, the use must be supervised by the guardian from the side, and the use of this product be able to guide and take responsibility.
8. This product is not suitable for cosmetic surgery, nose augmentation, etc. - cut the parts that have undergone skin beauty surgery, pregnant women, it is not recommended for patients with heart disease.