


LOMI
 CORDLESS
 WATER FLOSSER



User Guide
LOMB9955

* Please read instruction for use before use



INDEX

IMPORTANT SAFEGUARDS	1
GENERAL INFORMATION AND SAFETY INSTRUCTIONS	1
DESCRIPTION OF MODEL DENTAL SPA	2
BEFORE USE	2
USE OF THE DENTAL SPA	2
BASIC INFORMATION ABOUT USING OF THE IRRIGATOR	3
CLEAN OF THE IRRIGATOR	3
DEVICE SYMBOL	4
TECHNICAL SPECIFICATION	4
WARRANTY CERTIFICATION	4
COUPON OF WARRANTY	5

IMPORTANT SAFEGUARDS

When using electrical product, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING.


GENERAL INFORMATION AND SAFETY INSTRUCTIONS

1. Do not dismantle this product at any time other than when disposing of it.
2. Do not use outdoors or operate where aerosol product are being used or oxygen is being used or administered. This is a fire hazard.
3. Use of this product may interfere with the operation of pacemakers or other internal medical devices. If you use one of these devices you should consult your physician before using this irrigator.
4. Keep out of reach of children, this product should only be used by children supervised by adults. Keep jet nozzle attachments away from small children.
5. To avoid bacterial growth and contamination in water, always empty the water reservoir completely after use and clean the pump system by running the unit for a few seconds when the reservoir is emptied.
6. Never push the nozzle against teeth or gums.

DESCRIPTION OF MODEL DENTAL SPA

Dental spa uses a penetrating water jet to clean the spaces between your teeth and to reach hard to clean areas that conventional toothbrushes may struggle to reach. It will also massage the gums. Regular irrigation will assist in the prevention of gum disease and improve your dental and oral hygiene. Always use the irrigator after brushing your teeth. Please read the use instructions carefully before the first use.

1. On/Off and mode selection button (HI-LO)
2. Device
3. Irrigator jet tip
4. flip top water inlet

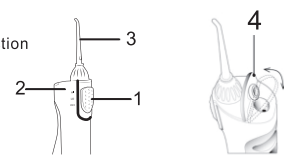

BEFORE USE

1. There may be moisture inside the main water unit on receipt of your irrigator. This is moisture left from the distilled water used to quality control the item before it is sold.
2. Insert 2 x AA batteries into the battery compartment at the base of the unit (see pic.). This compartment can be opened with a coin.
3. Close the battery cover securely after installing the batteries.



USE OF THE DENTAL SPA

1. Attach the irrigator jet tip to the device unit it locks in place (pic. 1).
2. Lift up the flip top on the device (pic. 2) and fill the reservoir with lukewarm water and add mouthwash, if desired (pic. 3) the water reservoir can be filled with cold or warm water (no hotter than 40°C), close the flip top firmly when the water reservoir beneath is full.
3. Insert the irrigator jet into your mouth, while you are leaning over a sink.
4. Hold the unit in an upright position with the irrigator jet tip pointing toward your teeth.
5. Gently push the switch upwards until it locks into the "Low" or "High" position. The pump will start operating (pic. 4). We recommend that you start off with the "Low" setting and move to the "High" setting as your gums get used to the massage motion and strengthen.




HIGH: For a thorough cleansing
LOW: For a delicate gums massage

6. Close lips enough to prevent splashing but keep your mouth slightly open to allow the water to run into the sink (pic. 5).



7. Move the jet along your teeth. Direct the jet stream at a 45 degree angle just above the gum line. Starting with the back teeth, follow the gum line and pause briefly before teeth and teeth. Continue until all areas around and between teeth have been cleaned (pic. 6).

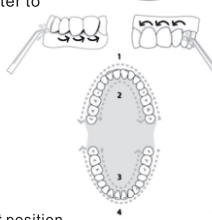

8. Ensure that you keep the device in an upright position. If you turn the irrigator upside down, it will not take water from the reservoir.
9. With the irrigator jet tip still attached, open the water reservoir by the top cap and drain away any remaining water/liquid.
10. Switch on the unit for a few seconds to pump out any water that remains in the pump reservoir or irrigator jet tip. We recommend that you do this after every use to avoid encouraging bacteria. Otherwise bacteria could collect in the remaining water. If the unit is not used for some time.
11. Should your dentist have recommended use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use with chlorhexidine, to prevent any deposit built-up.

BASIC INFORMATION ABOUT USING OF THE IRRIGATOR

1. Never fill the reservoir with hot water. Gums can be very sensitive and easily damaged by hot water.
2. Never turn the unit on before putting the irrigator jet in your mouth.
3. Never push the tip irrigator jet tip against teeth or gums.
4. Never use this device for cleaning any other part of the body such as ear or eyes.
5. Please note that those with sensitive gums may experience light bleeding the first few times they use this irrigator. If bleeding persists, please consult your doctor.
6. Place handle back on base and completely empty the tank after use.

CLEAN OF THE IRRIGATOR

1. Switch off the unit (1), remove the water reservoir by sliding it downward (2) and remove the irrigator jet tip.



2. Wash the reservoir in warm water (A). Wipe the white housing and the device with a damp cloth (B).



3. Never use caustic or abrasive cleaning agents as they may damage the device. Particular care should be taken when washing the hose connections in the reservoir.
4. Never submerge the unit in the water. After washing wipe with a dry cloth.

DEVICE SYMBOL



Attention



Correct disposal of the product (electric and electronic waste) (Applicable in the countries of the European Union and in the countries using systems for separate collection of waste). The symbol on the product or its manuals indicates that the product complies with the directive WEEE (waste electrical and electronic equipment) and must not be disposed of among domestic waste. At the end of its lifecycle, the user shall dispose of the equipment in specific collection areas. Otherwise, the sanctions foreseen by the current laws about waste are applicable. For detailed information concerning collection systems, refer to the local waste disposal service.

 According to EEC Directive: 2004/108/CE Electromagnetic Compatibility 2006/95/CE Directive Low Voltage

 Read carefully the instructions

TECHNICAL SPECIFICATION

Power source : 2x AA alkaline batteries

90-Day Limited Warranty

Territory: United States/Canada

LIMITED WARRANTY TO ORIGINAL CONSUMER

This product as supplied and distributed new by an authorized retailer is warranted by Southern Telecom, Inc. to the original consumer purchaser against defects in materials and workmanship ("Warranty") as follows:

- To obtain warranty service:**
- Visit our customer support portal: www.customersupport123.com
 - Select the brand and model of your device
 - Select "Service Request" and fill out the form to begin your inquiry
 - Exchanges require dated proof of purchase from an authorized retailer

Your product will be repaired or replaced, at our option, for the same or similar model of equal value if examination by the service center determines this product is defective. Products received damaged as a result of shipping will require you to file a claim with the carrier.

Warranty service not provided

This warranty does not cover damage resulting from accident, misuse, abuse, improper installation or operation, lack of reasonable care, unauthorized modification including software modifications such as the installation of custom firmware. This warranty is voided in the event any unauthorized person opens, alters or repairs this product. All products being returned to the authorized service center for repair must be suitably packaged.

NO WARRANTIES, WHETHER EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OTHER THAN THOSE EXPRESSLY DESCRIBED ABOVE SHALL APPLY. DISTRIBUTOR FURTHER DISCLAIMS ALL WARRANTIES AFTER THE EXPRESS WARRANTY PERIOD STATED ABOVE. NO OTHER EXPRESS WARRANTY OR GUARANTEE GIVEN BY ANY PERSON, FIRM OR ENTITY WITH RESPECT TO THE PRODUCT SHALL BE BINDING ON DISTRIBUTOR. REPAIR, REPLACEMENT, OR REFUND OF THE ORIGINAL PURCHASE PRICE - AT DISTRIBUTOR'S SOLE DISCRETION - ARE THE EXCLUSIVE REMEDIES OF THE CONSUMER. IN NO EVENT WILL DISTRIBUTOR OR ITS MANUFACTURERS BE LIABLE FOR ANY INCIDENTAL, DIRECT, INDIRECT, SPECIAL, PUNITIVE OR CONSEQUENTIAL DAMAGES (SUCH AS, BUT NOT LIMITED TO, DAMAGES FOR LOSS OF PROFITS, BUSINESS, SAVINGS, DATA OR RECORDS) CAUSED BY THE USE, MISUSE OR INABILITY TO USE THE PRODUCT. EXCEPT AS STATED HEREIN, NO OTHER WARRANTIES SHALL APPLY. NOTWITHSTANDING THE FOREGOING, CONSUMER'S RECOVERY AGAINST DISTRIBUTOR SHALL NOT EXCEED THE PURCHASE PRICE OF THE PRODUCT SOLD BY DISTRIBUTOR. THIS WARRANTY SHALL NOT EXTEND TO ANYONE OTHER THAN THE ORIGINAL CONSUMER WHO PURCHASED THE PRODUCT AND IS NOT TRANSFERABLE.

Some countries, states or provinces do not allow the exclusion or limitation of incidental or consequential damages or allow limitations on warranties, so limitation or exclusions may not apply to you. This warranty gives you specific legal rights, and you may have other rights that vary from state to state or province to province. Contact your authorized retailer to determine if another warranty applies.