


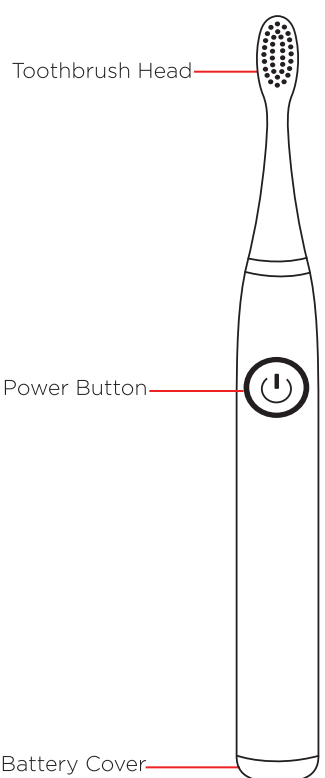





Lomi User Guide

(Front)

<p>Using the Toothbrush</p> <ol style="list-style-type: none"> 1. Install batteries, install a toothbrush head, and power ON. 2. Wet the brush head, and apply toothpaste. 3. Brush your teeth for at least 3 minutes, as recommended by dentists. If any pain or bleeding persists from brushing your teeth, contact your dentist. 4. When finished brushing, rinse the brush head only, then dry. <p>Changing the Brush Head Change the brush head after 3 months, or when the brush appears worn.</p> <p>Cleaning Wipe the outside of the toothbrush with a damp cloth or sponge. Do NOT submerge in water or other liquid.</p> <p>Storage For long-term storage, remove the battery and store in a cool, dry place away from sunlight, heat, wind and humidity.</p> <p>The toothbrush is best stored in an airtight box or bag (<i>not included</i>).</p> 	<p style="text-align: center;">Warranty Information</p> <p>90-Day Limited Warranty Territory: United States/Canada</p> <p>LIMITED WARRANTY TO ORIGINAL CONSUMER This product as supplied and distributed new by an authorized retailer is warranted by Southern Telecom, Inc. to the original consumer purchaser against defects in materials and workmanship ("Warranty") as follows:</p> <p>To obtain warranty service:</p> <ul style="list-style-type: none"> • Visit our customer support portal: www.customersupport123.com • Select the brand and model of your device • Select "Service Request" and fill out the form to begin your inquiry • Exchanges require dated proof of purchase from an authorized retailer <p>Your product will be repaired or replaced, at our option, for the same or similar model of equal value if examination by the service center determines this product is defective. Products received damaged as a result of shipping will require you to file a claim with the carrier.</p> <p>Warranty service not provided This warranty does not cover damage resulting from accident, misuse, abuse, improper installation or operation, lack of reasonable care, unauthorized modification including software modifications such as the installation of custom firmware. This warranty is voided in the event any unauthorized person opens, alters or repairs this product. All products being returned to the authorized service center for repair must be suitably packaged.</p> <p>NO WARRANTIES, WHETHER EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OTHER THAN THOSE EXPRESSLY DESCRIBED ABOVE SHALL APPLY. DISTRIBUTOR FURTHER DISCLAIMS ALL WARRANTIES AFTER THE EXPRESS WARRANTY PERIOD STATED ABOVE. NO OTHER EXPRESS WARRANTY OR GUARANTEE GIVEN BY ANY PERSON, FIRM OR ENTITY WITH RESPECT TO THE PRODUCT SHALL BE BINDING ON DISTRIBUTOR, REPAIR, REPLACEMENT, OR REFUND OF THE ORIGINAL PURCHASE PRICE - AT DISTRIBUTOR'S SOLE DISCRETION - ARE THE EXCLUSIVE REMEDIES OF THE CONSUMER. IN NO EVENT WILL DISTRIBUTOR OR ITS MANUFACTURERS BE LIABLE FOR ANY INCIDENTAL, DIRECT, INDIRECT, SPECIAL, PUNITIVE OR CONSEQUENTIAL DAMAGES (SUCH AS, BUT NOT LIMITED TO, DAMAGES FOR LOSS OF PROFITS, BUSINESS, SAVINGS, DATA OR RECORDS) CAUSED BY THE USE, MISUSE OR INABILITY TO USE THE PRODUCT, EXCEPT AS STATED HEREIN. NO OTHER WARRANTIES SHALL APPLY, NOTWITHSTANDING THE FOREGOING, CONSUMER'S RECOVERY AGAINST DISTRIBUTOR SHALL NOT EXCEED THE PURCHASE PRICE OF THE PRODUCT SOLD BY DISTRIBUTOR, THIS WARRANTY SHALL NOT EXTEND TO ANYONE OTHER THAN THE ORIGINAL CONSUMER WHO PURCHASED THE PRODUCT AND IS NOT TRANSFERABLE.</p> <p>Some countries, states or provinces do not allow the exclusion or limitation of incidental or consequential damages or allow limitations on warranties, so limitation or exclusions may not apply to you. This warranty gives you specific legal rights, and you may have other rights that vary from state to state or province to province. Contact your authorized retailer to determine if another warranty applies.</p> <p>All product names, trademarks, and registered trademarks are property of their respective owners. All company, product, and service names used in this document are for identification purposes only. Use of these names, trademarks, and brands does not imply endorsement.</p> <p style="text-align: right;">Printed in China Designed in the U.S. v1.06/2021</p>	  <p style="text-align: center;">User Guide LOMC1001</p>
---	--	---

<p style="text-align: center;">Before You Begin...</p> <p>READ CAREFULLY AND KEEP THIS MANUAL</p> <p>Caution</p> <ul style="list-style-type: none"> • THIS IS NOT A TOY, please keep away from children • Not recommended for use by children • Keep plastic bags and small parts away from children • Use this appliance only for its intended use as described in this manual • DO NOT use attachments not recommended by the manufacturer • To avoid injury, inspect the toothbrush for damage before using. Test parts of toothbrush gently with your hands. • Never operate this appliance if it is not working properly or if it has been dropped or damaged • Keep away from rain and humidity • Do NOT immerse in water. This can cause shock or injury • Keep away from open flames • Keep the appliance away from heated surfaces • Supervision is recommended when this appliance is used by, on or near individuals with certain disabilities • Do not operate after taking pain killers, sedatives or alcohol • Do not open this product or attempt to repair the unit yourself. Doing so could cause shock or injury, and will invalidate the manufacturer's warranty <p>Questions & Support If you have any questions or problems concerning this product, please visit our website at customersupport123.com.</p> <p>In the Box</p> <ul style="list-style-type: none"> • Electric Toothbrush (x1) • Brush Heads (x2) • AA Battery (x1) • Instruction Manual (x1) 	<p style="text-align: center;">Location of Controls</p> 	<p style="text-align: center;">Getting Started</p> <p>Installing Batteries</p> <ol style="list-style-type: none"> 1. Turn the battery cover on the bottom of the toothbrush counterclockwise to loosen, and remove. 2. Install the included AA battery into the battery compartment following the polarity (+/-) diagram in the compartment. 3. Replace the battery cover and turn clockwise until tight to close.  <p>Installing a Brush Head To install: Place the brush head on the top of the toothbrush, and push down all the way. Turn clockwise until the bristles are aligned with the power button to lock.</p>  <p>To detach: Turn the toothbrush head counterclockwise until it clicks, and gently pull the toothbrush head off.</p> <p>Powering ON/OFF</p> <ul style="list-style-type: none"> • Press the button on the front of the toothbrush to power ON or OFF.  <p>Instructions continued on the next page.</p>
--	--	---

(Back)