

MÄHLI

AQUA FLOSS

MANUAL

WATER FLOSSER



User Manual

MDB-1101

Before You Begin...

READ CAREFULLY AND KEEP THIS MANUAL

Caution

- **THIS IS NOT A TOY, please keep away from children**
- **Not recommended for use by children**
- Keep plastic bags and small parts away from children
- Use this appliance only for its intended use as described in this manual
- DO NOT use attachments not recommended by the manufacturer
- To avoid injury, inspect the aqua floss for damage before using. Test parts of aqua floss gently with your hands.
- Never operate this appliance if it is not working properly or if it has been dropped or damaged
- Keep away from open flames
- Keep the appliance away from heated surfaces
- Supervision is recommended when this appliance is used by, on or near individuals with certain disabilities
- Do not operate after taking pain killers, sedatives or alcohol
- Do not open this product or attempt to repair the unit yourself. Doing so could cause shock or injury, and will invalidate the manufacturer's warranty
- Please consult the local authorities on the correct way to dispose of this device.

Questions & Support

If you have any questions or problems concerning this product, please visit our website at customersupport123.com.


In the Box

- Water Flosser (x1)
- Nozzle
- User Manual (x1)

Specifications:

- Material: food grade ABS, food grade silicone
- Capacity: 100ml
- Net Weight: 110g

Location of Controls

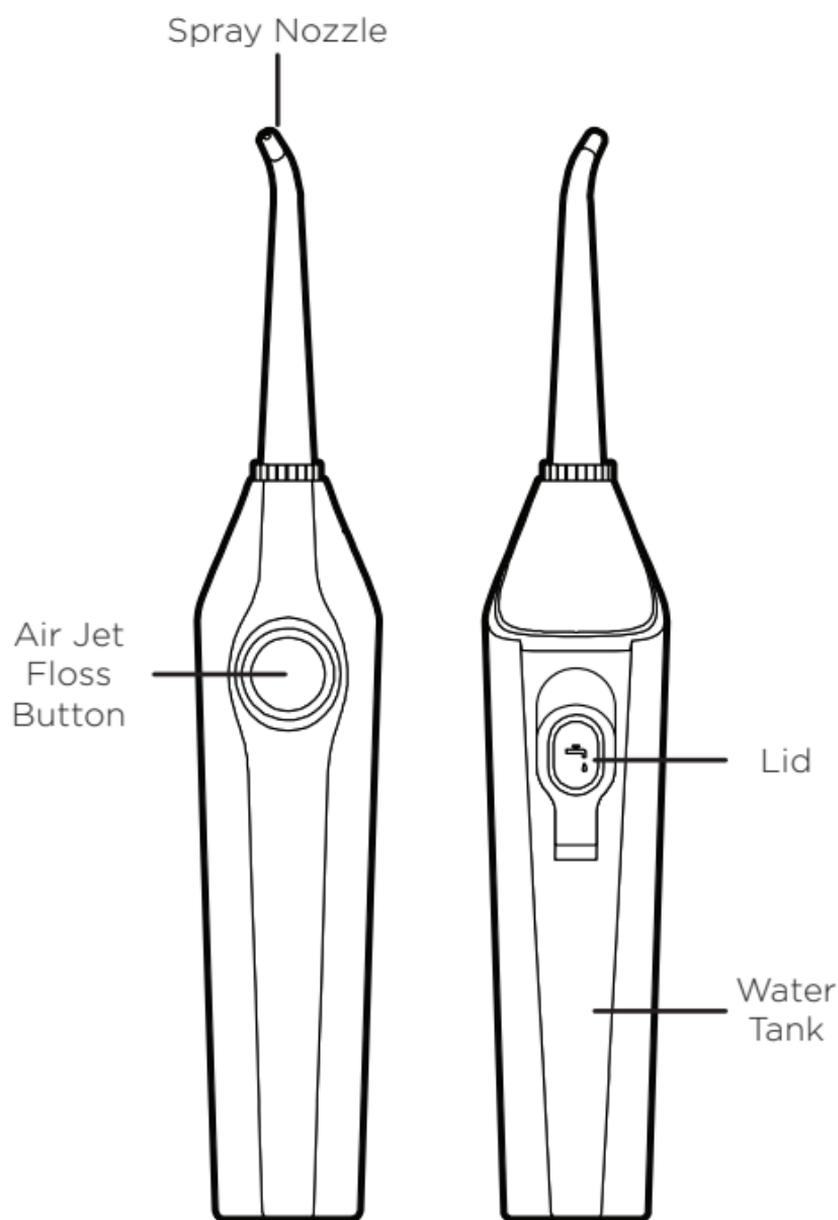

FRONT

BACK

Getting Started

Introduction:

Pressurized water stream cleans hard-to-reach areas in-between teeth and removes plaque that brushes leave behind, helping maintain oral hygiene for a healthier smile.

No batteries, no hassle. Easy-to-use one button design provides the most efficient manual water flossing method, simply fill the tank with water and press down to start flossing.

Highly portable design fits into a purse or carry-on, providing a perfect travel solution to maintaining oral cleanliness.

Quiet operation lets you clean your teeth without disturbing others, perfect for a cleaner smile any time of day or night.

How to use:

1. Open the lid, and pour the water into the water tank. The maximum capacity of the water tank is 100mL. Do not over fill past this level.




2. Press the air jet floss button to rinse the mouth.



3. Use it over the sink, keep lips slightly open to prevent water from splashing to allow the water to drip into the sink.



4. Keep the nozzle at about 45 degree angle along the teeth/gums, and move it along as the jet is activated. This will clean food debris between teeth and gums. Do this in short bursts flushing until all the teeth and gums are cleaned.

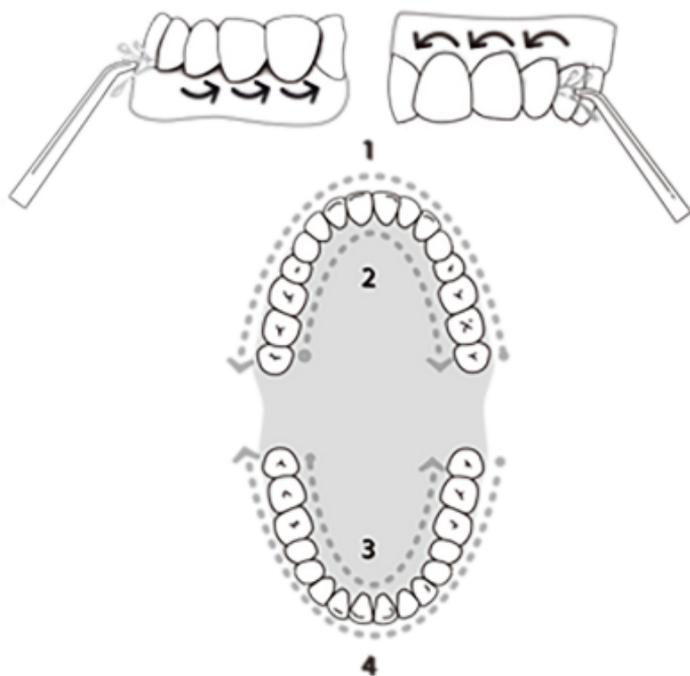

Note: This product is for oral cavity use only, NOT for eyes, ears or nose.



Rating Label Definitions



This symbol signifies that the product must not be discarded as household waste, and should be delivered to an appropriate collection facility for recycling. Proper disposal and recycling helps protect natural resources, human health and the environment.

Warranty Information

90-day Limited Warranty

Territory: United States/Canada

LIMITED WARRANTY TO ORIGINAL CONSUMER

This product (including any accessories included in the original packaging) as supplied and distributed new by an authorized retailer is warranted by Southern Telecom, Inc. to the original consumer purchaser against defects in materials and workmanship ("Warranty") as follows:

To obtain warranty service:

- Visit our customer support portal: www.customersupport123.com
- Select the brand and model of your device
- Select **Service Request** and fill out the form to begin your inquiry
- Exchanges require dated proof of purchase from an authorized retailer

Your product will be repaired or replaced, at our option, for the same or similar model of equal value if examination by the service center determines this product is defective. Products received damaged as a result of shipping will require you to file a claim with the carrier.

Warranty service not provided

This warranty does not cover damage resulting from accident, misuse, abuse, improper installation or operation, lack of reasonable care, unauthorized modification including software modifications such as the installation of custom firmware. This warranty is voided in the event any unauthorized person opens, alters or repairs this product. All products being returned to the authorized service center for repair must be suitably packaged.

NO WARRANTIES, WHETHER EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OTHER THAN THOSE EXPRESSLY DESCRIBED ABOVE SHALL APPLY. DISTRIBUTOR FURTHER DISCLAIMS ALL WARRANTIES AFTER THE EXPRESS WARRANTY PERIOD STATED ABOVE. NO OTHER EXPRESS WARRANTY OR GUARANTY GIVEN BY ANY PERSON, FIRM OR ENTITY WITH RESPECT TO THE PRODUCT SHALL BE BINDING ON DISTRIBUTOR. REPAIR, REPLACEMENT, OR REFUND OF THE ORIGINAL PURCHASE PRICE - AT DISTRIBUTOR'S SOLE DISCRETION - ARE THE EXCLUSIVE REMEDIES OF THE CONSUMER. IN NO EVENT WILL DISTRIBUTOR OR ITS MANUFACTURERS BE LIABLE FOR ANY INCIDENTAL, DIRECT, INDIRECT, SPECIAL, PUNITIVE OR CONSEQUENTIAL DAMAGES (SUCH AS, BUT NOT LIMITED TO, DAMAGES FOR LOSS OF PROFITS, BUSINESS, SAVINGS, DATA OR RECORDS) CAUSED BY THE USE, MISUSE OR INABILITY TO USE THE PRODUCT. EXCEPT AS STATED HEREIN, NO OTHER WARRANTIES SHALL APPLY. NOTWITHSTANDING THE FOREGOING, CONSUMER'S RECOVERY AGAINST DISTRIBUTOR SHALL NOT EXCEED THE PURCHASE PRICE OF THE PRODUCT SOLD BY DISTRIBUTOR. THIS WARRANTY SHALL NOT EXTEND TO ANYONE OTHER THAN THE ORIGINAL CONSUMER WHO PURCHASED THE PRODUCT AND IS NOT TRANSFERABLE.

Some countries, states or provinces do not allow the exclusion or limitation of incidental or consequential damages or allow limitations on warranties, so limitation or exclusions may not apply to you. This warranty gives you specific legal rights, and you may have other rights that vary from state to state or province to province. Contact your authorized retailer to determine if another warranty applies.

MÄHLI

All product names, trademarks, and registered trademarks are property of their respective owners. All company, product, and service names used in this document are for identification purposes only. Use of these names, trademarks, and brands does not imply endorsement.