

# MÄHLI

## **RECHARGEABLE SONIC PLUS** ELECTRIC TOOTHBRUSH



User Manual  
MDB-1115

# Before You Begin...

## READ CAREFULLY AND KEEP THIS MANUAL

### Caution

To reduce any risk of accident, take the following precautions:

- Never use the appliance while you are in the bath.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision.
- Do not manipulate the appliance.
- Regular care and cleaning is essential in maintaining your cleaning brush in peak working condition.
- Always clean the unit thoroughly after each use.
- DO NOT immerse in water or other liquid (water resistant IPX7).
- DO NOT use harsh cleaning products containing alcohol, petrol, or acetone.
- DO NOT use if you have a history of skin conditions or problems without first seeking medical advice.
- If you experience discomfort, discontinue use immediately and consult a physician.
- DO NOT use on sensitive or damaged skin.
- Close supervision is necessary when this device is used by/near children or disabled.
- Please keep these instructions in a safe place for future reference.
- Use this device only for its intended use as described in this manual.
- This product is suitable for cleaning teeth, gums, tongue and skin. Do not use it for other purpose.
- IMPORTANT: Avoid using in conjunction with abrasives to whiten and remove tooth stains.
- Check that the power transformer has the proper voltage.
- IMPORTANT: Do not use high resistance products on teeth and gums.
- Do not use in small children, those who cannot operate on their own or those with a weak oral cavity.
- IMPORTANT: If many of your teeth are being treated or the implant is not firm, do not use it.
- Don't drop it, don't hit it.
- Don't soak it in water for a long time.


### Questions & Support

If you have any questions or problems concerning this product, please visit our website at [customersupport123.com](http://customersupport123.com).

### In the Box

- Electric Toothbrush (x1)
- Brush Heads (x4)
- USB Charging Cable (x1)
- User Manual (x1)

## Location of Controls

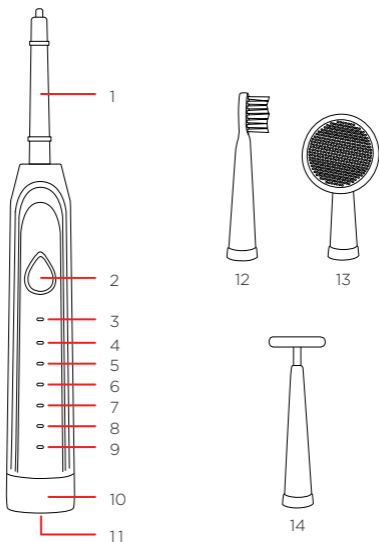

1. Interchangeable Head
2. Power ON/OFF Switch
3. Whitening Mode
4. Massage Mode
5. Cleaning Mode
6. Sensitive Teeth & Gums Mode
7. Polish Mode
8. Gum Care Mode

9. LED Charging Indicator
10. Waterproof Cover
11. USB Charging Port
12. Dental Care Brush Head
13. Silicone Brush Head
14. Facelift Bar Head

# Getting Started

## Charging the Toothbrush

Turn counterclockwise and remove the waterproof cover at the bottom of the toothbrush, and connect the small connector of the included USB charging cable into the USB charging port.

Connect the larger connect to a USB power source, such as the USB port on a computer, or a **5V USB wall adapter** (*not included*).

Once the toothbrush is fully charged, remove the USB charging cable from the toothbrush. And put the waterproof cover back on the bottom and turn it clockwise to tightly close it.


## Using the Toothbrush

This product has 6 working modes. Choose the desired working mode that best suits you.

- **Whitening Mode:** Whitening teeth and restoring beauty
- **Massage Mode:** Frequency conversion impact massage gum
- **Cleaning Mode:** Cleaning mouth with sound wave vibration
- **Soft Mode:** Gentle care for sensitive toothbrush
- **Polish Mode:** Polished and polished teeth
- **Gum Protection Mode:** Clean the four corners, clean the mouth

## Attaching a brush

Place the head at a 45° angle and turn it clockwise to align. To remove the head, turn it counterclockwise and then pull it up.

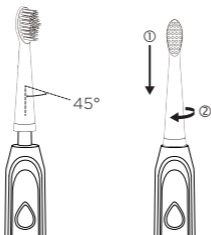

## Dental Cleaning

1. Moisten the brush and apply a suitable amount of toothpaste.
2. Place the brush on the teeth at a 45° angle to the gum line.
3. Press the power switch to start brushing your teeth.



45° angle for the gumline.

**Note:** Gently move the brush to the part to be cleaned, and the pressure should be within 150g (very low). The toothbrushes use the high-frequency vibration to clean the oral cavity, so the pressure required when brushing your teeth is smaller than that of normal toothbrushes. About 1/8 of the pressure of normal brushes is sufficient.

The dentist recommends that the oral cavity be divided into four areas. During the brushing process, move the toothbrush slowly to make sure each area has a cleaning time of approximately 30 seconds. This machine does not have a 15-second area reminder function and will stop for 0.5 every 15 seconds. To remind you to change areas to brush your teeth.

## Facial Cleansing

1. Lather your face and use the **Silicone Brush Head** to thoroughly cleanse your skin.
2. Press the power switch to turn it on and choose your suitable working mode.
3. Use the cleansing brush to trace small circles all over the face, emphasizing the oiliest areas that tend to produce more blackheads. Three minutes will be enough.



Silicone  
Brush Head

## Facelift

1. Apply your facial cream or facial treatment on your face, and attach the **Facelift Bar Head** to the toothbrush.
2. Press the power switch to turn it on and choose your suitable working mode.



Facelift  
Bar Head

## Cleaning

Wipe the outside of the toothbrush with a damp cloth or sponge. Do NOT submerge in water or other liquid.

## Storing

Store in a cool, dry place away from sunlight, heat, wind and humidity.

# Compliance Information

## FCC Information

This device complies with Part 15 of the FCC Rules.

**Caution:** Any changes or modifications not expressly approved could void the user's authority to operate the equipment.

**Note:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

## Radiation Exposure Information

This device has been evaluated to meet general RF exposure requirements. The device can be used in portable exposure conditions without restriction.

## WARNING Li-ion Battery Inside

This product is fitted with a Li-ion battery. Do not damage, open, or dismantle the battery and do not use it in damp and/or corrosive conditions. Use only with compatible chargers. Never dispose of batteries in a fire, and never expose them to high temperatures. Do not expose the product to temperatures exceeding 60°C (140°F).



## Rating Label Definitions



Direct current symbol.



This symbol signifies that the product must not be discarded as household waste, and should be delivered to an appropriate collection facility for recycling. Proper disposal and recycling helps protect natural resources, human health and the environment.



This symbol signifies and declares that this equipment complies with battery charging which complies with California CEC and energy efficiency.



This symbol represents FCC Declaration of Conformity and compliance.

All product names, trademarks, and registered trademarks are property of their respective owners. All company, product, and service names used in this document are for identification purposes only. Use of these names, trademarks, and brands does not imply endorsement.

Printed in China  
Designed in the U.S.  
v1. 06/2021

# Warranty Information

## 90-day Limited Warranty

Territory: United States/Canada

### LIMITED WARRANTY TO ORIGINAL CONSUMER

This product (including any accessories included in the original packaging) as supplied and distributed new by an authorized retailer is warranted by Southern Telecom, Inc. to the original consumer purchaser against defects in materials and workmanship ("Warranty") as follows:

#### To obtain warranty service:

- Visit our customer support portal: [www.customersupport123.com](http://www.customersupport123.com)
- Select the brand and model of your device
- Select **Service Request** and fill out the form to begin your inquiry
- Exchanges require dated proof of purchase from an authorized retailer

Your product will be repaired or replaced, at our option, for the same or similar model of equal value if examination by the service center determines this product is defective. Products received damaged as a result of shipping will require you to file a claim with the carrier.

#### Warranty service not provided

This warranty does not cover damage resulting from accident, misuse, abuse, improper installation or operation, lack of reasonable care, unauthorized modification including software modifications such as the installation of custom firmware. This warranty is voided in the event any unauthorized person opens, alters or repairs this product. All products being returned to the authorized service center for repair must be suitably packaged.

NO WARRANTIES, WHETHER EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OTHER THAN THOSE EXPRESSLY DESCRIBED ABOVE SHALL APPLY. DISTRIBUTOR FURTHER DISCLAIMS ALL WARRANTIES AFTER THE EXPRESS WARRANTY PERIOD STATED ABOVE. NO OTHER EXPRESS WARRANTY OR GUARANTY GIVEN BY ANY PERSON, FIRM OR ENTITY WITH RESPECT TO THE PRODUCT SHALL BE BINDING ON DISTRIBUTOR. REPAIR, REPLACEMENT, OR REFUND OF THE ORIGINAL PURCHASE PRICE - AT DISTRIBUTOR'S SOLE DISCRETION - ARE THE EXCLUSIVE REMEDIES OF THE CONSUMER. IN NO EVENT WILL DISTRIBUTOR OR ITS MANUFACTURERS BE LIABLE FOR ANY INCIDENTAL, DIRECT, INDIRECT, SPECIAL, PUNITIVE OR CONSEQUENTIAL DAMAGES (SUCH AS, BUT NOT LIMITED TO, DAMAGES FOR LOSS OF PROFITS, BUSINESS, SAVINGS, DATA OR RECORDS) CAUSED BY THE USE, MISUSE OR INABILITY TO USE THE PRODUCT. EXCEPT AS STATED HEREIN, NO OTHER WARRANTIES SHALL APPLY. NOTWITHSTANDING THE FOREGOING, CONSUMER'S RECOVERY AGAINST DISTRIBUTOR SHALL NOT EXCEED THE PURCHASE PRICE OF THE PRODUCT SOLD BY DISTRIBUTOR. THIS WARRANTY SHALL NOT EXTEND TO ANYONE OTHER THAN THE ORIGINAL CONSUMER WHO PURCHASED THE PRODUCT AND IS NOT TRANSFERABLE.

Some countries, states or provinces do not allow the exclusion or limitation of incidental or consequential damages or allow limitations on warranties, so limitation or exclusions may not apply to you. This warranty gives you specific legal rights, and you may have other rights that vary from state to state or province to province. Contact your authorized retailer to determine if another warranty applies.