

MÄHLI

CORDLESS HANDHELD MASSAGER

**HANDHELD DEEP TISSUE
& MUSCLE MASSAGER**



User Manual
MDM-55

Before You Begin...

READ CAREFULLY AND KEEP THIS MANUAL

Caution

- **THIS IS NOT A TOY, please keep away from children**
- **Not recommended for use by children**
- Keep plastic bags and small parts away from children
- Use this appliance only for its intended use as described in this manual
- DO NOT use attachments not recommended by the manufacturer
- To avoid injury, inspect the toothbrush for damage before using. Test parts of toothbrush gently with your hands.
- Never operate this appliance if it is not working properly or if it has been dropped or damaged
- Keep away from rain and humidity
- Do NOT immerse in water. This can cause shock or injury
- Keep away from open flames
- Keep the appliance away from heated surfaces
- Supervision is recommended when this appliance is used by, on or near individuals with certain disabilities
- Do not operate after taking pain killers, sedatives or alcohol
- Do not open this product or attempt to repair the unit yourself. Doing so could cause shock or injury, and will invalidate the manufacturer's warranty
- This device contains a built-in battery and should not be disposed of with regular household waste. Please consult the local authorities on the correct way to dispose of this device.


Questions & Support

If you have any questions or problems concerning this product, please visit our website at customersupport123.com.

In the Box

- Massage Gun (x1)
- Vibration Heads (x3)
- USB Charging Cable (x1)
- User Manual (x1)

Location Of Controls

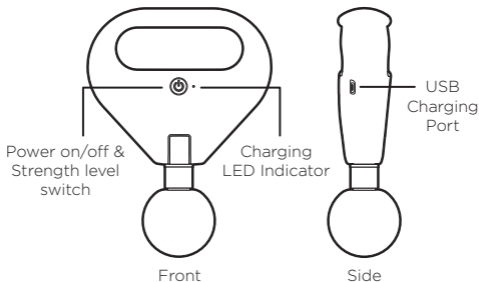

Removable, replaceable vibration heads, 3 heads in total.
*When removing the vibration head, just exert your strength to pull it out.



Ball head
(For large muscle groups)




Small round head
(For joints)



U head
(For vertebrae)

Charging

1. Plug the small connector of the supplied USB charging cable to the USB charging port. Plug the large connector into the USB port on a computer or into a compatible **5V USB adapter** (not included) and plug into an electrical outlet.

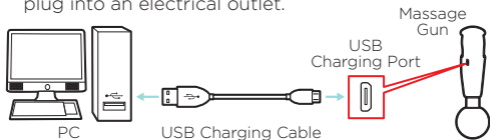

Getting Started

2. The charging LED indicator in the adapter will turn RED while charging and turn OFF when fully charged.

Note: Fully charge the massage gun before using it for the first time. It is not recommended to run out the battery before charging.

Changing the vibration Heads

Firmly hold the massager and pull out the vibration head from the massager. Attach another vibration head to the socket of the massager.

Using the Massager

Choose a vibration head that best suits the area you wish to massage. Attach the head securely to the massager.

- Press and hold the Power button to power on.
- Press the Power button to adjust the speeds. There are two speeds in total.

Apply the massager head to your body, and move the massager around the area in a circular motion.

Caution: In order to reduce the risk of the electric shock, fire and personal injury, this product must be used in accordance with the following instructions.

- For adults only
- It can only be used on dry, clean skin surfaces of the body, and cannot be used through clothes. Press and move it gently on the skin. The duration on each part is about 60 seconds.
- Use this device only on soft tissues of the body without causing pain or discomfort. Do not use it on head or any hard or boned part of the body.

Getting Started

- Only use the vibrating head that can produce the best results.
- Bruise may be caused whether under moderate pressure or normal pressure. Please always check the area that receives physiotherapy.
- Stop using it immediately once there is any sign of pain or discomfort.
- Keep fingers, hair or other parts of the body away from the shaft and back of the vibrating head, because they might be pinched.
- Do not immerse it in water or allow water.
- Do not drop the device or use it improperly,
- Carefully check the device and battery before each use.
- Do not charge or modify the device in any way.
- Never operate or charge the device in an unattended way.

Compliance Information

FCC Information

This device complies with Part 15 of the FCC Rules.

Caution: Any changes or modifications not expressly approved could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Radiation Exposure Information

This device has been evaluated to meet general RF exposure requirements. The device can be used in portable exposure conditions without restriction.

WARNING Li-ion Battery Inside

This product is fitted with a Li-ion battery. Do not damage, open, or dismantle the battery and do not use it in damp and/or corrosive conditions. Use only with compatible chargers. Never dispose of batteries in a fire, and never expose them to high temperatures. Do not expose the product to temperatures exceeding 60°C (140°F).



Rating Label Definitions



Direct current symbol.



This symbol signifies that the product must not be discarded as household waste, and should be delivered to an appropriate collection facility for recycling. Proper disposal and recycling helps protect natural resources, human health and the environment.



This symbol signifies and declares that this equipment complies with battery charging which complies with California CEC and energy efficiency.



This symbol represents FCC Declaration of Conformity and compliance.

Warranty Information

90-day Limited Warranty

Territory: United States/Canada

LIMITED WARRANTY TO ORIGINAL CONSUMER

This product (including any accessories included in the original packaging) as supplied and distributed new by an authorized retailer is warranted by Southern Telecom, Inc. to the original consumer purchaser against defects in materials and workmanship ("Warranty") as follows:

To obtain warranty service:

- Visit our customer support portal: www.customersupport123.com
- Select the brand and model of your device
- Select **Service Request** and fill out the form to begin your inquiry
- Exchanges require dated proof of purchase from an authorized retailer

Your product will be repaired or replaced, at our option, for the same or similar model of equal value if examination by the service center determines this product is defective. Products received damaged as a result of shipping will require you to file a claim with the carrier.

Warranty service not provided

This warranty does not cover damage resulting from accident, misuse, abuse, improper installation or operation, lack of reasonable care, unauthorized modification including software modifications such as the installation of custom firmware. This warranty is voided in the event any unauthorized person opens, alters or repairs this product. All products being returned to the authorized service center for repair must be suitably packaged.

NO WARRANTIES, WHETHER EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, OTHER THAN THOSE EXPRESSLY DESCRIBED ABOVE SHALL APPLY. DISTRIBUTOR FURTHER DISCLAIMS ALL WARRANTIES AFTER THE EXPRESS WARRANTY PERIOD STATED ABOVE. NO OTHER EXPRESS WARRANTY OR GUARANTY GIVEN BY ANY PERSON, FIRM OR ENTITY WITH RESPECT TO THE PRODUCT SHALL BE BINDING ON DISTRIBUTOR. REPAIR, REPLACEMENT, OR REFUND OF THE ORIGINAL PURCHASE PRICE - AT DISTRIBUTOR'S SOLE DISCRETION - ARE THE EXCLUSIVE REMEDIES OF THE CONSUMER. IN NO EVENT WILL DISTRIBUTOR OR ITS MANUFACTURERS BE LIABLE FOR ANY INCIDENTAL, DIRECT, INDIRECT, SPECIAL, PUNITIVE OR CONSEQUENTIAL DAMAGES (SUCH AS, BUT NOT LIMITED TO, DAMAGES FOR LOSS OF PROFITS, BUSINESS, SAVINGS, DATA OR RECORDS) CAUSED BY THE USE, MISUSE OR INABILITY TO USE THE PRODUCT. EXCEPT AS STATED HEREIN, NO OTHER WARRANTIES SHALL APPLY. NOTWITHSTANDING THE FOREGOING, CONSUMER'S RECOVERY AGAINST DISTRIBUTOR SHALL NOT EXCEED THE PURCHASE PRICE OF THE PRODUCT SOLD BY DISTRIBUTOR. THIS WARRANTY SHALL NOT EXTEND TO ANYONE OTHER THAN THE ORIGINAL CONSUMER WHO PURCHASED THE PRODUCT AND IS NOT TRANSFERABLE.

Some countries, states or provinces do not allow the exclusion or limitation of incidental or consequential damages or allow limitations on warranties, so limitation or exclusions may not apply to you. This warranty gives you specific legal rights, and you may have other rights that vary from state to state or province to province. Contact your authorized retailer to determine if another warranty applies.

MÄHLI

All product names, trademarks, and registered trademarks are property of their respective owners. All company, product, and service names used in this document are for identification purposes only. Use of these names, trademarks, and brands does not imply endorsement.

Printed in China
Designed in the U.S.
v1. 06/2021